



Feel the NEW Standard of Care^{*}

The PowerPort[®] Implanted Port is an implantable access device designed to provide repeated access to the vascular system. Port access is performed by percutaneous needle insertion using a non-coring needle. Power injection is performed using a PowerLoc[®] Safety Infusion Set only.



The PowerPort[®] system can be used for infusion of medications, I.V. fluids, parenteral nutrition solutions, blood products, and for the withdrawal of blood samples. All materials are biocompatible, can be used with virtually all injectable solutions and are safe with CECT.


What's Different About The PowerPort[®] Implanted Port?

The PowerPort[®] Implanted Port is the **FIRST** implantable access device indicated for power injection of contrast media during CECT scans.

- Contrast enhanced computed tomography (CECT) produces quick, accurate visualizations used to track tumor markers or pulmonary embolisms.
- During power injection, the CT Team injects a contrast agent at a high-rate into the bloodstream in order to achieve the most detailed images of the area being scanned.
- Power injection produces superior images that can reveal small details in soft tissue.
- Like all implanted ports, the PowerPort[®] device eliminates repeated needle sticks in the arm or wrist.

IDENTIFYING A PATIENT WITH THE POWERPORT[®] IMPLANTED PORT


- Check patient's chart for a PowerPort[®] device patient **record sticker**.
- Palpate top of port to identify **three palpation points** (bumps) on the septum, arranged in a triangle.
- Palpate the sides of the port to identify **triangular** port housing.
- **Ask the patient.** Patients with a PowerPort[®] Implanted Port should have a patient identification card, ID bracelet or key ring to help remind them they have a PowerPort[®] device.
- **Always verify** the patient has a PowerPort[®] device by at least two means, and ensure it is accessed with a PowerLoc[®] Safety Infusion Set, prior to power injection.

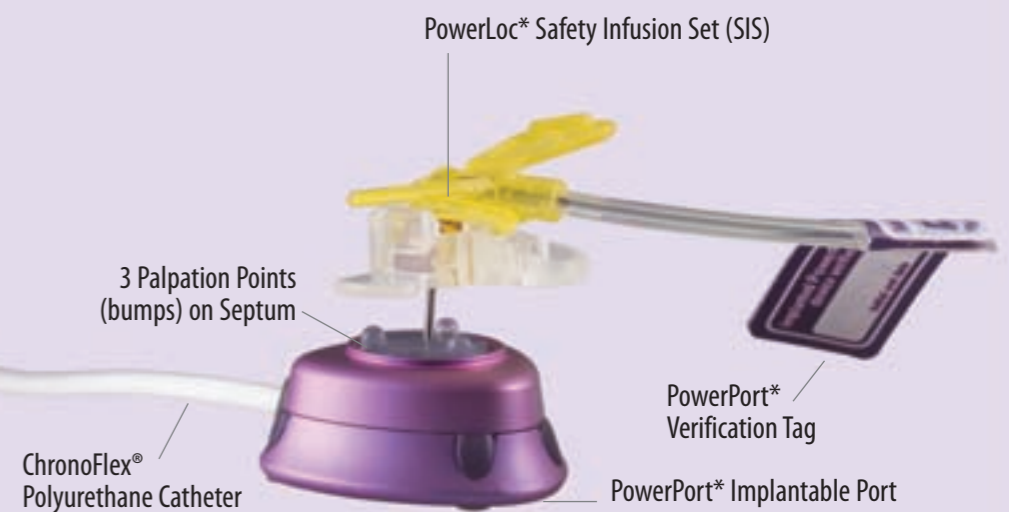

Check patient chart



Feel for bumps on septum

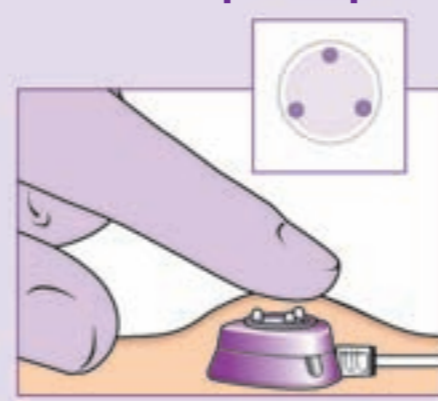

Feel for triangle shape



Ask your patient



Double Check




USE AND MAINTENANCE


Accessing Implanted Ports

Procedure:

1. Perform aseptic site preparation.
2. Locate port septum by palpation.
 - a. Locate base of port with non-dominant hand.
 - b. Triangulate port between thumb and first two fingers of non-dominant hand. Aim for center point of these three fingers.

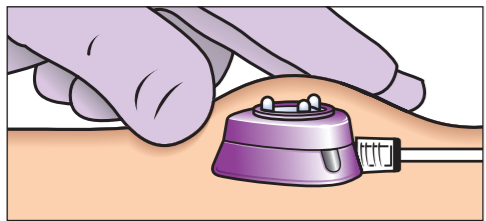

3. Insert PowerLoc[®] Safety Infusion Set perpendicular to port septum. Advance needle through the skin and septum until reaching bottom of reservoir.

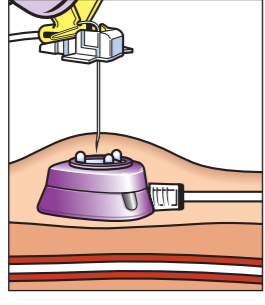

4. Confirm correct positioning of the needle within the port reservoir by aspiration of blood ("flashback"). If there is doubt regarding proper needle placement, have a radiographic dye procedure done to confirm placement.
5. Always flush the port following injection.
6. Perform heparin lock procedure for open-ended catheters. **Caution:** Remember that some patients may be hypersensitive to heparin or suffer from heparin induced thrombocytopenia (HIT) and these patients must not have their port locked with heparinized saline.
7. After therapy completion, flush port per institutional protocol. Close clamp while injecting the last 0.5 ml of flush solution. Use positive pressure technique.

Heparin Lock

To help prevent clot formation and catheter blockage, implanted ports with open-ended catheters should be flushed with 10 ml sterile normal saline using a turbulent push-pause flushing method

after each use followed by 5 ml of heparinized saline. Clamp the tubing while infusing the last 0.5 ml of fluid to reduce potential for blood back-flow into the catheter tip, which could encourage catheter clotting. If the port remains unused for long periods of time, the 5 ml heparin solution should be changed at least every four weeks. **Caution:** Remember that some patients may be hypersensitive to heparin or suffer from heparin induced thrombocytopenia (HIT) and these patients must not have their port locked with heparinized saline.

Determine Port Volume

For PowerPort[®] implanted ports, you will need to check the patient's chart to determine the length of catheter used for each individual patient. For PowerPort[®] catheters, multiply the catheter length in cm by 0.02 ml, then add 0.60 ml for the port reservoir:

Example:
Catheter length: _____ cm x 0.02 ml/cm
+0.60 ml (port septum) = _____ ml volume,
total priming volume for patient port and catheter.

Recommended Flushing Volumes:

Procedure	Volume (100 U/ml)
When port not in use	5 ml heparinized saline every 4 weeks
After each infusion of medication or TPN	10 ml sterile normal saline, then 5 ml heparinized saline
After blood withdrawal	20 ml sterile normal saline, then 5 ml heparinized saline
After power injection of contrast media	10 ml sterile normal saline, then 5 ml heparinized saline

Procedure:

Review Site Preparation in the PowerPort[®] Nursing Guide, and Accessing Implanted Ports section before proceeding with the following:

1. Explain procedure to patient and prepare injection site.
2. Attach a syringe filled with sterile normal saline to needle.

3. Aseptically locate and access port with PowerLoc[®] Infusion Set, or other non-coring safety needle. Confirm correct positioning of the needle within the port reservoir by aspiration of blood ("flashback"). If there is doubt regarding proper needle placement, have a radiographic dye procedure done to confirm placement.

4. After therapy completion, flush port per institutional protocol, then repeat with 5 ml 100 U/ml heparinized saline, or with volume calculated above. Close clamp while injecting last 0.5 ml of flush solution.
5. To deaccess PowerLoc[®] safety infusion set from the port, activate safety mechanism while withdrawing needle until you feel a "click" at which time the needle should be captured within the safety mechanism of the PowerLoc[®] safety infusion set.

Note: Alcohol should not be used to soak or decontaminate polyurethane catheters because alcohol is known to degrade the polyurethane catheters over time with repeated and prolonged exposure.

SEE POWERPORT[®] NURSING GUIDE FOR MORE DETAILS.

CHECK FOR PATENCY PRIOR TO POWER INJECTION

Aspirate for adequate blood return and vigorously flush the port with at least 10 ml of sterile normal saline.

WARNING:

Failure to ensure patency of the catheter prior to power injection studies may result in port system failure.

TROUBLESHOOTING GUIDE

I. Aspiration Difficulties: DO NOT POWER INJECT IF YOU CANNOT ASPIRATE AS PATIENT INJURY MAY RESULT

A. Possible Causes:

1. Failure to flush adequately, resulting in lumen obstruction.
2. Catheter tip sucking up to vein wall with aspiration.
3. Blood clot, fibrin sheath, or particulate matter obstructing lumen when catheter is aspirated.
 - A clot or other obstruction in the catheter lumen can produce a one-way valve effect. During infusion, the catheter wall expands slightly and allows fluid to flow around the plug. During aspiration, the catheter wall contracts slightly, tightening down around the obstruction and preventing aspiration.
 - Fibrin sheaths usually begin to form within a few days after the insertion of a central venous catheter. If it has grown enough to extend the tip of the catheter, it may be pulled into and obstruct the catheter opening when aspiration is attempted, but no resistance to infusion.
4. Compression or transection of the catheter between the clavicle and first rib ("pinch-off area").
5. Kinked catheter.
 - Catheter may be pulled too tightly through skin tunnel, causing kink at vessel insertion site, or where it curves into the subcutaneous tunnel.
 - Catheter may be curled or kinked within the vessel, or under the dressing.
6. Malposition of catheter tip (i.e. jugular vein, outside of vein).
7. Improper catheter length selection for patient size.

B. Possible Solutions:

1. If no resistance to infusion is felt, attempt to flush with 10 ml normal saline. Then pull back gently on syringe plunger 2-3 ml, pause and proceed with aspiration.
2. If resistance to infusion is felt, check for signs of extravasation. If present, notify physician of possible catheter leakage or transection and embolization. If not present, see step 4.
3. Attempt to aspirate with a 20 ml syringe.
4. Move patient's arm, shoulder and head to see if a change in position will allow aspiration. If aspiration can only be accomplished with the patient in a certain position, the patient should be examined to see if the catheter has been placed in the "pinch-off" area.
5. Obtain physician's order for a chest x-ray to determine the position of the catheter.
 - If the catheter tip is not in the superior vena cava, the catheter should be repositioned.
 - If the catheter tip is not in a vein, the catheter should be replaced.
 - If the catheter has been placed through the "pinch-off" area, between the clavicle and the first rib, and is being compressed enough to interfere with infusion or aspiration, it is at risk for catheter transection and embolization. The physician should evaluate the patient for catheter replacement.

II. Patient with Fever Infection:

Symptoms:

- Inflammation at incision site
 - Fever
 - Positive site culture / or blood cultures
- If signs of infection are present:**
- Notify physician

III. Insufficient Flow: DO NOT POWER INJECT IF RESISTANCE TO FLUSHING SEEMS EXCESSIVE

Excessive force must not be used to flush an obstructed lumen. Insufficient blood flow may be caused by the catheter contacting the wall of the vein. The physician may attempt to dissolve the clot with a fibrinolytic agent before power injecting. Physician discretion advised.

Equipment:

- PowerLoc[®] Safety Infusion Set, or other non-coring safety needle. Choose a needle length based on reservoir depth, tissue thickness and the thickness of any dressing beneath the bend of the needle.
- Syringe containing port priming volume of a fibrinolytic agent.
- Syringe filled with sterile normal saline.

Procedure:

Review Site Preparation in the PowerPort[®] Nursing Guide and Accessing Implanted Ports section before proceeding with the following:

1. Explain procedure to patient and prepare injection site.
2. Aseptically locate and access the desired septum with needle attached to syringe, void of air and filled with port priming volume of fibrinolytic agent.

Warning:

If accessing a PowerPort[®] device with a PowerLoc[®] Safety Infusion Set, do not affix the PowerLoc[®] sticker that indicates the system can be power injected. Power injecting a blocked catheter could lead to catheter damage and patient injury.

3. Gently instill fibrinolytic solution. Use a gentle pull-push action on the syringe plunger to maximize solution mixing within port and catheter. **Warning:** Occluded catheter may not accept all of the solution. If strong resistance is felt, do not attempt to force into catheter.
4. Leave solution in place according to drug manufacturers recommendation and/or doctor's orders.
5. Attempt to aspirate solution and the clot(s).
6. If the clot(s) cannot be aspirated, repeat procedure.
7. Once the blockage has been aspirated and discarded, flush catheter with at least 20 ml of sterile normal saline.
8. Flush the catheter with 5 ml of heparinized saline. Remember that some patients may be hypersensitive to heparin, or suffer from heparin induced thrombocytopenia (HIT). These patients must not have their ports flushed with heparinized saline.
9. After therapy completion, flush port per institutional protocol. Close clamp while injecting last 0.5 ml of flush solution. Use positive pressure technique.

IV. Catheter Occlusion: DO NOT POWER INJECT AN OCCLUDED DEVICE

A. Possible Causes:

1. Blood clot completely obstructing lumen.
2. May be kinked, coiled, damaged, or compressed between the clavicle and the first rib.
3. Catheter tip may not be within vein.

V. Use of Fibrinolytic Agent for Catheter Blockage

Use of a fibrinolytic agent has successfully cleared clotted catheters when gentle irrigation and aspiration have failed. The instructions provided by the drug manufacturer should be followed. Alcohol should not be used to soak or decontaminate polyurethane catheters because alcohol is known to degrade polyurethane catheters over time with repeated and prolonged exposure.

Signs of Pinch-off:

- Clinical:**
- Difficulty with blood withdrawal
 - Resistance to infusion of fluids
 - Patient position changes required for infusion of fluids or blood withdrawal
- Radiologic:**
- Grade 1 or 2 distortion on chest X-ray.
- "Pinch-off" should be evaluated for degree of severity prior to explanation. Patients indicating any degree of catheter distortion at the clavicle/first rib area should be followed diligently. There are grades of "pinch-off" that should be recognized with appropriate chest x-ray as follows:

Grade	Severity	Recommended Action
Grade 0	No distortion	No action
Grade 1	Distortion present without luminal narrowing	Chest x-ray should be taken every one to three months to monitor progression of pinch off to grade 2 distortion. Shoulder position during chest x-rays should be noted as it can contribute to changes in distortion grades.
Grade 2	Distortion present with luminal narrowing	Removal of the catheter should be considered.
Grade 3	Catheter transection or fracture	Prompt removal of the catheter.

Please consult product labels and inserts for any indications, contraindications, hazards, warnings, cautions and instructions for use.